

# TMT10plex Mass Tag Labeling Kits and Reagents

90110 90111 90406 90113

2457.2

Number	Description
90110	<b>TMT10plex™ Isobaric Label Reagent Set</b> , sufficient reagents for one 10plex isobaric experiment <b>Contents:</b> TMT <sup>10</sup> -126™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -127N™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -127C™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -128N™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -128C™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -129N™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -129C™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -130N™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -130C™ Label Reagent, 1 × 0.8mg TMT <sup>10</sup> -131™ Label Reagent, 1 × 0.8mg
90111	<b>TMT10plex Isobaric Label Reagent Set</b> , sufficient reagents for three 10plex isobaric experiments <b>Contents:</b> TMT <sup>10</sup> -126 Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -127N Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -127C Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -128N Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -128C Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -129N Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -129C Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -130N Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -130C Label Reagent, 3 × 0.8mg TMT <sup>10</sup> -131 Label Reagent, 3 × 0.8mg
90406	<b>TMT10plex Isobaric Label Reagent Set</b> , sufficient reagents for one 10plex isobaric experiment <b>Contents:</b> TMT <sup>10</sup> -126 Label Reagent, 1 × 5mg TMT <sup>10</sup> -127N Label Reagent, 1 × 5mg TMT <sup>10</sup> -127C Label Reagent, 1 × 5mg TMT <sup>10</sup> -128N Label Reagent, 1 × 5mg TMT <sup>10</sup> -128C Label Reagent, 1 × 5mg TMT <sup>10</sup> -129N Label Reagent, 1 × 5mg TMT <sup>10</sup> -129C Label Reagent, 1 × 5mg TMT <sup>10</sup> -130N Label Reagent, 1 × 5mg TMT <sup>10</sup> -130C Label Reagent, 1 × 5mg TMT <sup>10</sup> -131 Label Reagent, 1 × 5mg

**90113 TMT10plex Isobaric Mass Tag Labeling Kit**, sufficient reagents for three 10plex isobaric experiments

**Contents:**

**TMT<sup>10</sup>-126 Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-127N Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-127C Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-128N Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-128C Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-129N Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-129C Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-130N Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-130C Label Reagent**, 3 × 0.8mg  
**TMT<sup>10</sup>-131 Label Reagent**, 3 × 0.8mg  
**Dissolution Buffer** (1M triethyl ammonium bicarbonate), 5mL  
**Denaturing Reagent** (10% SDS), 1mL  
**Reducing Reagent** (0.5M TCEP), 1mL  
**Iodoacetamide**, 12 × 9mg  
**Quenching Reagent** (50% hydroxylamine), 1mL  
**Pierce Trypsin Protease, MS Grade**, 5 × 20µg  
**Trypsin Storage Solution**, 250µL  
**Albumin, Bovine**, 2.5mg

**Storage:** Upon receipt store at -20°C. Reagents are shipped with dry ice.

**Note:** Products are for research use only – do not use for diagnostic procedures.


## Table of Contents

Introduction .....	2
Procedure Summary.....	3
Important Product Information .....	4
Additional Materials Required.....	4
Material Preparation .....	5
Preparing and Labeling Peptides with the TMT Isobaric Mass Tags .....	5
Troubleshooting .....	6
Additional Information .....	6
A. Data Acquisition Methods .....	6
B. Data Analysis and Quantitation .....	7
C. Information Available from Our Website .....	7
Related Thermo Scientific Products .....	8
General References .....	8
Limited Use Label License: TMT10plex™ Isobaric Mass Tag Labeling Kits and Reagents.....	8

## Introduction

The Thermo Scientific™ TMT™ Mass Tag Labeling Kits and Reagents enable multiplex relative quantitation by mass spectrometry (MS). Each mass-tagging reagent within a set has the same nominal mass (i.e., isobaric) and chemical structure composed of an amine-reactive NHS-ester group, a spacer arm and a mass reporter (Figure 1). The reagent set can be used to label up to ten different peptide samples prepared from cells or tissues. For each sample, a unique reporter mass (i.e., TMT<sup>10</sup> 126-131Da) in the low mass region of the MS/MS spectrum is used to measure relative protein expression levels during peptide fragmentation.

The Thermo Scientific™ TMT10plex™ Label Reagents share an identical structure with Thermo Scientific™ TMTzero™ and TMTsixplex™ Reagents but contain different numbers and combinations of  $^{13}\text{C}$  and  $^{15}\text{N}$  isotopes in the mass reporter. The different isotopes result in a 10plex set of tags that have monoisotopic mass differences in the reporter that can be detected using high resolution Thermo Scientific™ Orbitrap™ Mass Spectrometry Instruments. Advantages of the TMT10plex Label Reagents include increased sample multiplexing for relative quantitation, increased sample throughput and fewer missing quantitative channels among samples.

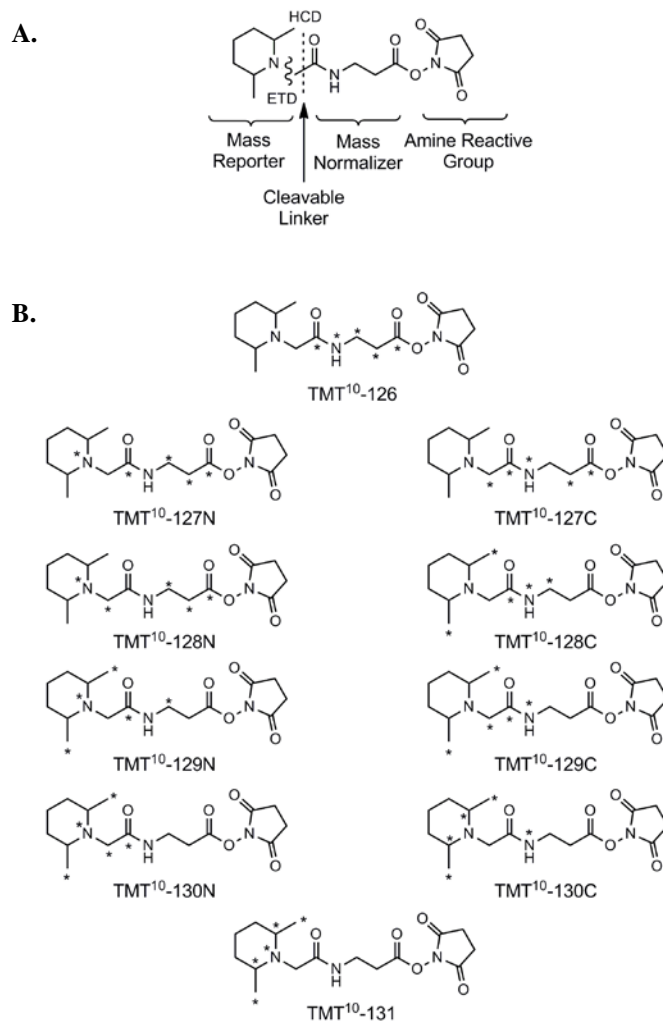


**Figure 1. Chemical structures of the Thermo Scientific TMT10plex Label Reagents.**

A. Functional regions of the reagent structure including MS/MS fragmentation sites by higher energy collision dissociation (HCD) and electron transfer dissociation (ETD). B. TMT10plex Reagent structures and isotope positions (\*).

## Procedure Summary

Protein extracts isolated from cells or tissues are reduced, alkylated and digested overnight. Samples are labeled with the TMT Reagents and then mixed before sample fractionation and clean-up. Labeled samples are analyzed by high resolution Orbitrap LC-MS/MS before data analysis to identify peptides and quantify reporter ion relative abundance (Figure 2).

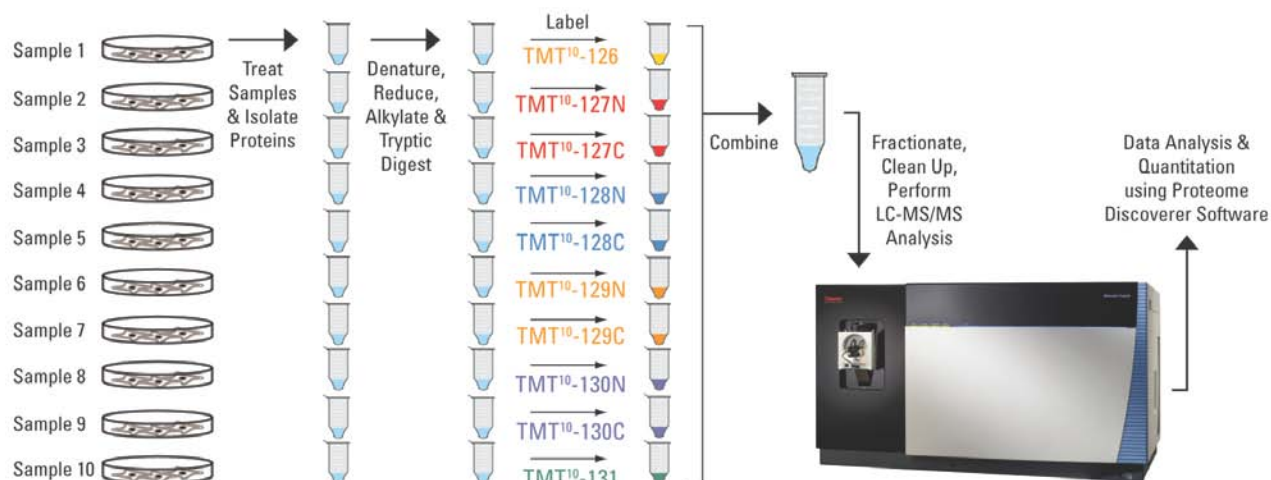

**Figure 2. Procedure schematic for using the Thermo Scientific TMT10plex Label Reagents.**

## Important Product Information

- The TMT Reagents are moisture-sensitive. To avoid moisture condensation onto the product, the vial must be equilibrated to room temperature before opening.
- The TMT Reagents are amine-reactive and modify lysine residues and peptide N-termini. All amine-containing buffers and additives must be removed before digestion and labeling.
- All samples must be digested, labeled and then mixed equally before desalting, fractionation and LC-MS/MS. For optimal results, use 25-100µg of peptide for each labeling reaction.
- To avoid contamination of MS samples, always wear gloves when handling samples and gels. Use ultrapure MS-grade reagents. Perform sample preparation in a clean work area.
- The TMTzero Label Reagent (Product No. 90067) can be used to optimize methods before multiplexed analysis of samples with TMT10plex Label Reagent sets.

## Additional Materials Required

- Microcentrifuge tubes
- Anhydrous acetonitrile (Acetonitrile, LC-MS Grade, Product No. 51101)
- Water, LC-MS Grade (Product No. 51140)
- Chilled (-20°C) acetone
- Protein assay (e.g., Thermo Scientific™ BCA Protein Assay Kit, Product No. 22235)
- 75-300µm capillary C<sub>18</sub> reversed-phase column
- High-resolution Orbitrap Mass Spectrometer with online liquid chromatography system (see Table 1 for recommended instruments)
- Data analysis software (e.g., Thermo Scientific™ Proteome Discoverer™ Software)
- Optional: C18 spin tips or columns (e.g., Thermo Scientific™ Pierce™ C18 Spin Columns, Product No. 89870 or Pierce™ C18 Tips, Product No. 87784)

## Material Preparation

**Note:** The 50% hydroxylamine and 10% SDS stock solutions provided with the kit may precipitate during storage. Warm both solutions to room temperature and vortex before use. The amounts listed below are sufficient for preparing and labeling 10 samples.

100mM TEAB (triethyl ammonium bicarbonate)	Add 500μL of the Dissolution Buffer (1M TEAB) to 4.5mL of ultrapure water.
Lysis Buffer	Add 200μL of the Denaturing Reagent (10% SDS) to 1.8mL of 100mM TEAB.
200mM TCEP	Add 70μL of the Reducing Reagent (0.5M TCEP) to 70μL of ultrapure water. Then add 35μL of the Dissolution Buffer (1M TEAB).
5% Hydroxylamine	Add 50μL of the Quenching Reagent (50% hydroxylamine) to 450μL of 100mM TEAB.

## Preparing and Labeling Peptides with the TMT Isobaric Mass Tags

**Note:** BSA can be used as a control sample for method optimization. Dissolve BSA to 1mg/mL using 100mM TEAB. Use 25-100μg of protein per labeling reaction. The Thermo Scientific™ Pierce™ Mass Spec Sample Prep Kit for Cultured Cells can also be used to prepare peptide digests for TMT reagent labeling.

### A. Preparing Whole Cell Protein Extracts

1. Culture cells to harvest at least 100μg of protein per condition. For best results, culture a minimum of  $2 \times 10^6$  cells.

**Note:** Rinse cells 2-3 times with 1X PBS to remove cell culture media. Pellet cells using low-speed centrifugation (i.e.,  $< 1000 \times g$ ) to prevent premature cell lysis.

2. Lyse the cells by adding five cell-pellet volumes of Lysis Buffer (i.e., 100μL of Lysis Buffer for a 20μL cell pellet).

**Note:** Lysis buffers such as Thermo Scientific™ RIPA Lysis and Extraction Buffer (Product No. 89901) or 8M urea (Product No. 29700) in 50mM TEAB or HEPES buffer, pH 8 may be used as alternative denaturing cell lysis buffers. For urea-based lysis buffer, protein samples must be diluted to  $< 1M$  urea before digestion, and the final C18 desalting step (C.6) is not optional. Addition of protease and/or phosphatase inhibitors during lysis is optional and may interfere with MS analysis.

**Note:** Depending on the Lysis Buffer used it may be necessary to reduce sample viscosity by shearing DNA using a microtip sonicator or addition of a nuclease (e.g., Thermo Scientific™ Pierce™ Universal Nuclease for Cell Lysis, Product No. 88700)

3. Centrifuge lysate at  $16,000 \times g$  for 10 minutes at 4°C.
4. Carefully separate the supernatant and transfer into a new tube.
5. Determine the protein concentration of the supernatant using established methods such as the BCA Protein Assay Kit (Product No. 23227).

**Note:** Use samples at  $\geq 2\text{mg/mL}$ . Less concentrated samples may be used; however, it might be necessary to use larger volumes of reducing/alkylating reagents.

6. Transfer 100μg per condition into a new tube and adjust to a final volume of 100μL with 100mM TEAB.
7. Add 5μL of the 200mM TCEP and incubate sample at 55°C for 1 hour.
8. Immediately before use, dissolve one tube of iodoacetamide (9mg) with 132μL of 100mM TEAB to make 375mM iodoacetamide. Protect solution from light.
9. Add 5μL of the 375mM iodoacetamide to the sample and incubate for 30 minutes protected from light at room temperature.
10. Add six volumes (~600μL) of pre-chilled (-20°C) acetone. Allow the precipitation to proceed overnight.
11. Centrifuge the samples at  $8000 \times g$  for 10 minutes at 4°C. Carefully invert the tubes to decant the acetone without disturbing the white pellet. Allow the pellet to dry for 2-3 minutes.

## B. Protein Digestion

1. Resuspend 100µg of acetone-precipitated (or lyophilized) protein pellets with 100µL of 100mM TEAB.  
**Note:** An acetone-precipitated pellet may not completely dissolve; however, after proteolysis at 37°C, all the protein (peptides) will be solubilized.
2. Immediately before use, add 20µL of the Trypsin Storage Solution to the bottom of the trypsin glass vial and incubate for 5 minutes. Store any remaining reagent in single-use volumes at -80°C (e.g., 2.5µg of trypsin per 100µg of protein).
3. Add 2.5µL of trypsin (i.e., 2.5µg) per 100µg of protein. Digest the sample overnight at 37°C.

## C. Peptide Labeling

1. Immediately before use, equilibrate the TMT Label Reagents to room temperature. For the 0.8mg vials, add 41µL of anhydrous acetonitrile to each tube. For the 5mg vials, add 256µL of solvent to each tube. Allow the reagent to dissolve for 5 minutes with occasional vortexing. Briefly centrifuge the tube to gather the solution.  
**Note:** Reagents dissolved in anhydrous acetonitrile are stable for one week when stored at -20°C and warmed to room temperature before opening. Anhydrous ethanol can be used as an alternative solvent to dissolve reagents.
2. Carefully add 41µL of the TMT Label Reagent to each 25-100µg sample. Alternatively, transfer the reduced and alkylated protein to the TMT Reagent vial.  
**Note:** A 100µL glass syringe or positive displacement pipette may be necessary to accurately measure and dispense TMT Reagents in volatile acetonitrile solvent.
3. Incubate the reaction for 1 hour at room temperature.
4. Add 8µL of 5% hydroxylamine to the sample and incubate for 15 minutes to quench the reaction.
5. Combine samples in a new microcentrifuge tube at equal amounts and store at -80°C.
6. Optional: Clean-up samples with C18 spin tips (Product No. 87784) or columns (Product No. 89870) before high-resolution LC-MS analysis. Peptide clean up is recommended before LC-MS analysis but is not required.

## Troubleshooting

Problem	Possible Cause	Solution
Poor labeling	A primary amine-based buffer was used (e.g., Tris, glycine)	Use non-primary amine-based buffers (e.g., TEAB, HEPES)
	Incorrect buffer pH	Make sure the buffer pH is ~8.0
	Too much sample was used	Label 25-100µg per sample
Protein precipitation	Lack of detergent present	Add detergent (e.g., 0.05-0.1% SDS) to the preparation
	pH decreased	Make sure pH is > 7.5

## Additional Information

### A. Data Acquisition Methods

Quantitation of peptides labeled with Thermo Scientific™ Tandem Mass Tag™ Reagents requires a high-resolution Orbitrap Mass Spectrometer capable of MS/MS fragmentation (Table 1). To resolve near-isobaric reporter ions, MS/MS resolution must be > 50,000 at 150 *m/z*. Higher energy collision dissociation (HCD) is recommended for TMT10plex reporter ion fragmentation. Optimal HCD fragmentation energy is instrument-dependent and can be optimized using TMTzero Reagents. Electron transfer dissociation (ETD) may be used as an alternative fragmentation method for peptide identification and quantitation; however, ETD is not recommended for TMT10plex Reagents because of reporter ion overlap (Table 2).

**Table 1. Instruments and MS/MS fragmentation options for peptide identification and quantitation with Thermo Scientific TMT Reagents.**

<u>Instrument</u>	<u>Fragmentation Method</u>	<u>Minimum Resolution Setting</u>	<u>Reference(s)</u>
Thermo Scientific Orbitrap Fusion™ Tribrid™ Mass Spectrometer	HCD/SPS-MS3	60,000	Viner, <i>et al.</i> (2013)
Thermo Scientific Orbitrap Elite™ Mass Spectrometer	HCD/MS3	30,000	Viner, <i>et al.</i> (2012)
Thermo Scientific Q Exactive™ Mass Spectrometer	HCD/MS2	35,000	Wühr, <i>et al.</i> (2012)
Thermo Scientific Orbitrap Velos Pro Mass Spectrometer	HCD/MS2	30,000	Ting, <i>et al.</i> (2011), Wenger, <i>et al.</i> (2011)

## B. Data Analysis and Quantitation

The peptide mass modification by the TMT10plex Reagents is identical to TMTsixplex Reagents and present in the UNIMOD database ([www.unimod.org](http://www.unimod.org)) and are listed below. Proteome Discoverer Software (1.4 and above) is recommended for TMT10plex relative quantitation. Additional software programs that may be used for TMT quantitation include Matrix Science™ Mascot™ Software (2.2 and above) and Proteome Software™ Scaffold™ Q+ Software. For data acquired using a combination of fragmentation methods (i.e., HCD/MS3 or HCD/ETD), Proteome Discoverer Software may be necessary to merge search results.

**Table 2. Modification masses of the Thermo Scientific TMT Label Reagents.**

<u>Label Reagent</u>	<u>Label Reagent</u>	<u>Modification Mass (monoisotopic)</u>	<u>Modification Mass (average)</u>	<u>HCD Monoisotopic Reporter Mass*</u>	<u>ETD Monoisotopic Reporter Mass**</u>
TMT <sup>10</sup> -126	TMT <sup>6</sup> -126	229.162932	229.2634	126.127725	114.127725
TMT <sup>10</sup> -127N	TMT <sup>6</sup> -127	229.162932	229.2634	127.124760	115.124760
TMT <sup>10</sup> -127C	-	229.162932	229.2634	127.131079	114.127725
TMT <sup>10</sup> -128N	-	229.162932	229.2634	128.128114	115.124760
TMT <sup>10</sup> -128C	TMT <sup>6</sup> -128	229.162932	229.2634	128.134433	116.134433
TMT <sup>10</sup> -129N	TMT <sup>6</sup> -129	229.162932	229.2634	129.131468	117.131468
TMT <sup>10</sup> -129C	-	229.162932	229.2634	129.137787	116.134433
TMT <sup>10</sup> -130N	-	229.162932	229.2634	130.134822	117.131468
TMT <sup>10</sup> -130C	TMT <sup>6</sup> -130	229.162932	229.2634	130.141141	118.141141
TMT <sup>10</sup> -131	TMT <sup>6</sup> -131	229.162932	229.2634	131.138176	119.138176

\* HCD is a collisional fragmentation method that generates ten unique reporter ions from 126 to 131Da.

\*\*ETD is a non-ergodic fragmentation method that generates six unique reporter ions from 114 to 119Da.

## C. Information Available from Our Website

- Tech Tip #49: Acetone precipitation of proteins
- Tech Tip #19: Remove detergent from protein samples



## Related Thermo Scientific Products

90114	1M Triethylammonium bicarbonate (TEAB), 50mL
90115	50% Hydroxylamine, 5mL
90067	TMTzero Label Reagent, 5 × 0.8mg
90061	TMTsixplex Isobaric Label Reagent Set, 1 × 0.8mg
90064	TMTsixplex Isobaric Mass Tagging Kit
90100	iodoTMTzero™ Label Reagent, 5 × 0.2mg
90101	iodoTMTsixplex™ Label Reagent Set, 1 × 0.2mg
90103	iodoTMTsixplex Isobaric Mass Tag Labeling Kit
84840	Pierce™ Mass Spec Sample Prep Kit for Cultured Cells
23227	BCA Protein Assay Kit
90057	Pierce Trypsin Protease, MS Grade
90051	Lys-C Protease, MS Grade
88300	Fe-NTA Phosphopeptide Enrichment Kit
88301	Pierce TiO <sub>2</sub> Phosphopeptide Enrichment and Clean-up Kit
89870	Pierce C18 Spin Columns, 25 columns
28904	Trifluoroacetic Acid, Sequanal Grade

## General References

- Altelaar A.F., *et al.* (2012). Benchmarking stable isotope labeling based quantitative proteomics. *J Proteomics* Oct 22. pii: S1874-3919(12)00704-X. doi: 10.1016/j.jprot.2012.10.009.
- Dillon, R., *et al.* (2011). Discovery of a novel B-Raf fusion protein related to c-Met drug resistance. *J Proteome Res* **10(11)**:5084-94.
- Lee, M.V., *et al.* (2011). A dynamic model of proteome changes reveals new roles for transcript alteration in yeast. *Mol Syst Biol* **7**: 514.
- McAllister, G.C., *et al.* (2012). Increasing the multiplexing capacity of TMTs using reporter ion isotopologues with isobaric masses. *Anal Chem* **84(17)**:7469-78.
- Ting, L., *et al.* (2011). MS3 eliminates ratio distortion in isobaric multiplexed quantitative proteomics. *Nat Methods* **8**: 937–40.
- Wenger, C.D., *et al.* (2011). Gas-phase purification enables accurate, multiplexed proteome quantification with isobaric tagging. *Nat Methods* **8(11)**:933-5.
- Werner, T., *et al.* (2012). High-resolution enabled TMT8-plexing. *Anal Chem* **84(16)**:7188-94.
- Wu, L., *et al.* (2013). Variation and genetic control of protein abundance in humans. *Nature* **499(7456)**:79-82.
- Wühr, M., *et al.* (2012). Accurate multiplexed proteomics at the MS2 level using the complement reporter ion cluster. *Anal Chem* **84(21)**:9214-21.
- Viner, R.I., *et al.* (2013). Increasing the multiplexing of protein quantitation from 6- to 10-Plex with reporter ion isotopologues. PN\_ASMS\_W617\_RViner\_R1.
- Viner, R.I., *et al.* (2012). Relative quantitation of TMT-labeled proteomes – Focus on sensitivity and precision. Application note #566.
- Viner, R.I., *et al.* (2009). Quantification of post-translationally modified peptides of bovine  $\alpha$ -crystallin using tandem mass tags and electron transfer dissociation. *J Proteomics* **72(5)**:874-85.
- Xiong, L., *et al.* (2011). Mass spectrometric studies on epigenetic interaction networks in cell differentiation. *J Biol Chem* **286(15)**:13657-68.
- Zhang, T., *et al.* (2010). Improving quantitation of TMT-labeled peptides using stepped higher-energy collisional dissociation. Application note #483.

## Limited Use Label License: TMT10plex™ Isobaric Mass Tag Labeling Kits and Reagents

(Thermo Scientific Product Nos. 90110, 90111, 90113 and 90406).

The purchase of this product conveys to the buyer the non-exclusive, non-transferable right to use the product and components of the product for internal research purposes of the buyer only, whether the buyer is an academic or for-profit entity. This product is not available for use in diagnostic procedures.

The buyer must obtain a license for use of the product in performing commercial services or for any other commercial purposes. For information on purchasing a license, contact proteomics.licensing@thermofisher.com, Pierce Biotechnology, Inc., 3747 North Meridian Road, Rockford, IL 61101.

If the buyer is not willing to accept the limitations of this limited use statement, Pierce Biotechnology, Inc. is willing to accept return of the unused product with a full refund. By opening and using this product, the buyer agrees to be bound by the conditions of this Limited Use Label License.



Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Product documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample.

NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED, INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS.

Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to humans or animals.

Current product instructions are available at [www.thermoscientific.com/pierce](http://www.thermoscientific.com/pierce). For a faxed copy, call 800-874-3723 or contact your local distributor.

© 2014 Thermo Fisher Scientific Inc. All rights reserved. Tandem Mass Tag, TMT, TMTzero, TMTsixplex, TMT10plex, iodoTMTzero and iodoTMTsixplex are trademarks of Proteome Sciences plc. Unless otherwise indicated, all other trademarks are property of Thermo Fisher Scientific Inc. and its subsidiaries. Printed in the USA.